



## New! Multi-Volume XL Pipette Tip


- 2-300  $\mu\text{L}$ —non-filtered volume range.
- 2-200  $\mu\text{L}$ —filtered volume range.

### One tip does it all!

- Graduations at 10, 50, 100, 150, and 200  $\mu\text{L}$ —eliminate the need to keep multiple tips in your inventory.
- Universal fit for broad compatibility with most major pipettor brands.
- Filter placement is optimized to accommodate 20  $\mu\text{L}$ , 100  $\mu\text{L}$ , and 200  $\mu\text{L}$  pipettors; non-filtered tips are also compatible with 300  $\mu\text{L}$  pipettors.

### Complete sample recovery!

- 91 mm (3.58 in) in length—easily reach the bottom of 12 x 75 mm tubes without risk of the pipettor nose cone touching the interior of the tube.
- Manufactured from virgin polypropylene to ensure purity and provide excellent optical clarity.
- NoStick<sup>®</sup> resin technology creates a low-binding surface for maximum sample recovery.
- Certified free of DNase, RNase, DNA, and PCR inhibitors.
- Sterile product is also pyrogen-free.
- Available in a standard 8 x 12 rack format.
- Patented press-to-open hinged lid enables easy, one-hand operation.

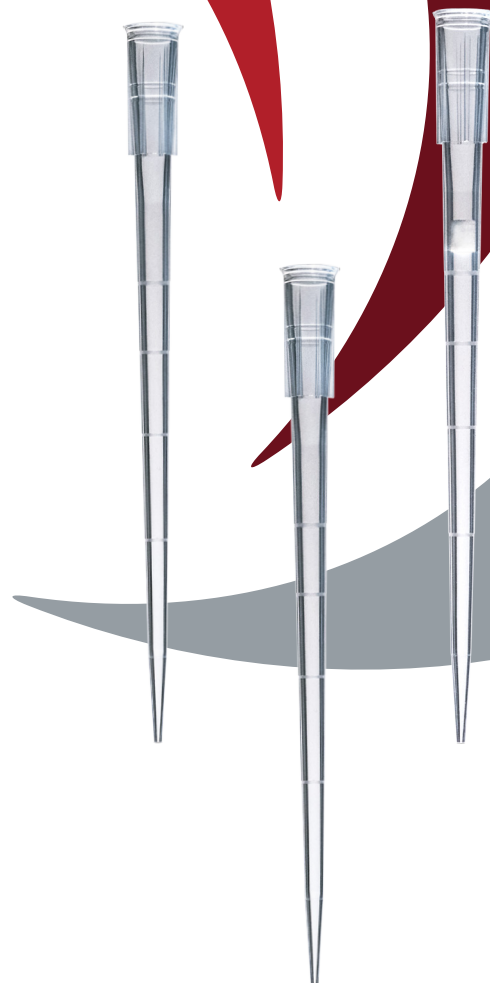

## SSIbio Multi-Volume XL Pipette Tip

Item #	Max Volume	Description	Color	Packaging
4280N00	300 $\mu\text{L}$	Multi-Volume XL Bulk	Natural	1000/pack, 10 packs/case
4287N00	300 $\mu\text{L}$	Multi-Volume XL Racked	Natural	2 racks of 96/unit, 5 units/case
4287NS0	300 $\mu\text{L}$	Multi-Volume XL Racked, Sterile	Natural	2 racks of 96/unit, 5 units/case
4287NSF	200 $\mu\text{L}$	Multi-Volume XL Filtered, Racked, Sterile	Natural	2 racks of 96/unit, 5 units/case